



Bolzano, 23rd February 2015

Subject: Informative regarding the modification of the indicator lights on Boscarol suction units

We are glad to inform our customers that the new display electronic circuits for the suction units OB 2012 and OB 1000 came into production.

The new indicator lights are intended to improve the intuitiveness of alerts, replacing the red colour, commonly associated to a situation of danger that might wrongly recall the idea of a malfunctioning.

- 1
- 2


OB 2012



The suction unit OB 2012 will have four yellow indicators (1) to signal the level of battery charge and a charge indicator (2), near the other ones, which signals with a yellow led the charging in progress and with a green led the completed charge.

OB 1000

The suction unit OB 1000 will have three indicator lights to signal (1): green for maximum autonomy, yellow for intermediate autonomy and orange for low autonomy or low battery; the charging indicator (2), next to the other leds, signals with a yellow led the charging in progress and with a green led the completed charge.

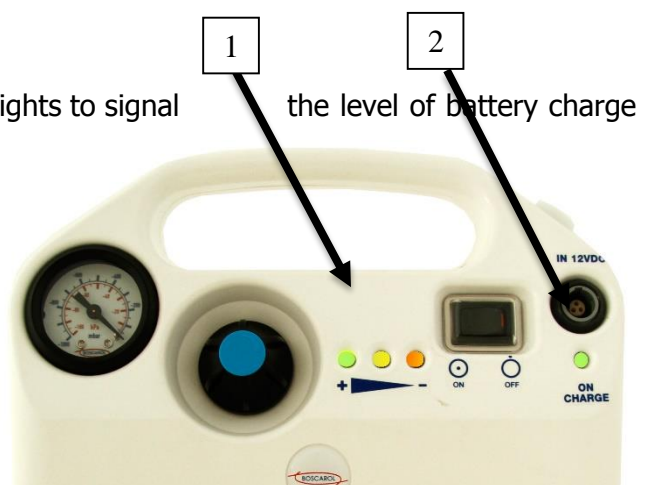



The light indicators will be placed, as previously, on the frontal part of the device, ensuring a good visibility also in case of direct light on the device.

The variation is only valid for the suction units that have been produced starting from august 2014, not on the devices already in use: it is in fact an improvement not explicitly requested by the notifying agency but, in our opinion, an initiative that will certainly make Boscarol suction units more user-friendly.

We remain at your disposal for any further assistance

Best regards

Oscar Boscarol srl